

Dragon® NaturallySpeaking® Medical 9.5 and Higher

Creating Macro Recorder Voice Commands

The Macro Recorder allows you to record a sequence of mouse movements and keystrokes. When you use one of the commands you created, the mouse movements and keystrokes you recorded are played back.

To Create a Macro Recorder Voice Command:


1. If it is not already running on your computer, start *Dragon NaturallySpeaking*.
2. On the **DragonBar**, click **Tools**.
3. Select **Add New Command** or say “Add new command.” The **MyCommands Editor** dialog box appears.
4. Make sure the cursor is in the **MyCommand Name** field.
5. Using your microphone, say a name for your command. For example, to insert a new contact into your EMR, you can create a command called “insert new contact.”
6. **(Optional)** In the **Description** field, dictate a short description of the command's function.
7. **(Optional but recommended)** From the **Group** dropdown list, select a group, or create a new group.
 - Note:* Although placing a command in a group is optional, doing so makes it easier for you to find commands in the **Command Browser** under the **Manage** mode.
 - Note:* To change a group name simply dictate or type over **User-defined** with a different group name. The new group name will appear in the list the next time you create a command.
8. For **Availability**, select **Application-specific**.
9. Open the application in which you want the command to work. For example, to have the Macro Recorder command work in your EMR, open the EMR application. Then return to the **MyCommands Editor** dialog box.
10. In the **Application** field, click the drop-down list to view all open applications.
11. Find your application and select it.



12. From the **Command Type** drop-down list, select **Macro Recorder**.
13. Click **Record** to begin recording. The **MyCommands Editor** window minimizes automatically. When the **Floating Recorder** window appears in the lower left corner of the screen, recording is already in process.
14. When you have finished recording, click the **Stop** button on the **Floating Recorder**.
15. Click the **X** in the upper right corner of the **Floating Recorder** to close the recorder and redisplay the **MyCommands Editor** dialog box.



Note: The recorded movements and/or keystrokes are shown in the **Actions** section of the **MyCommands Editor** dialog box.

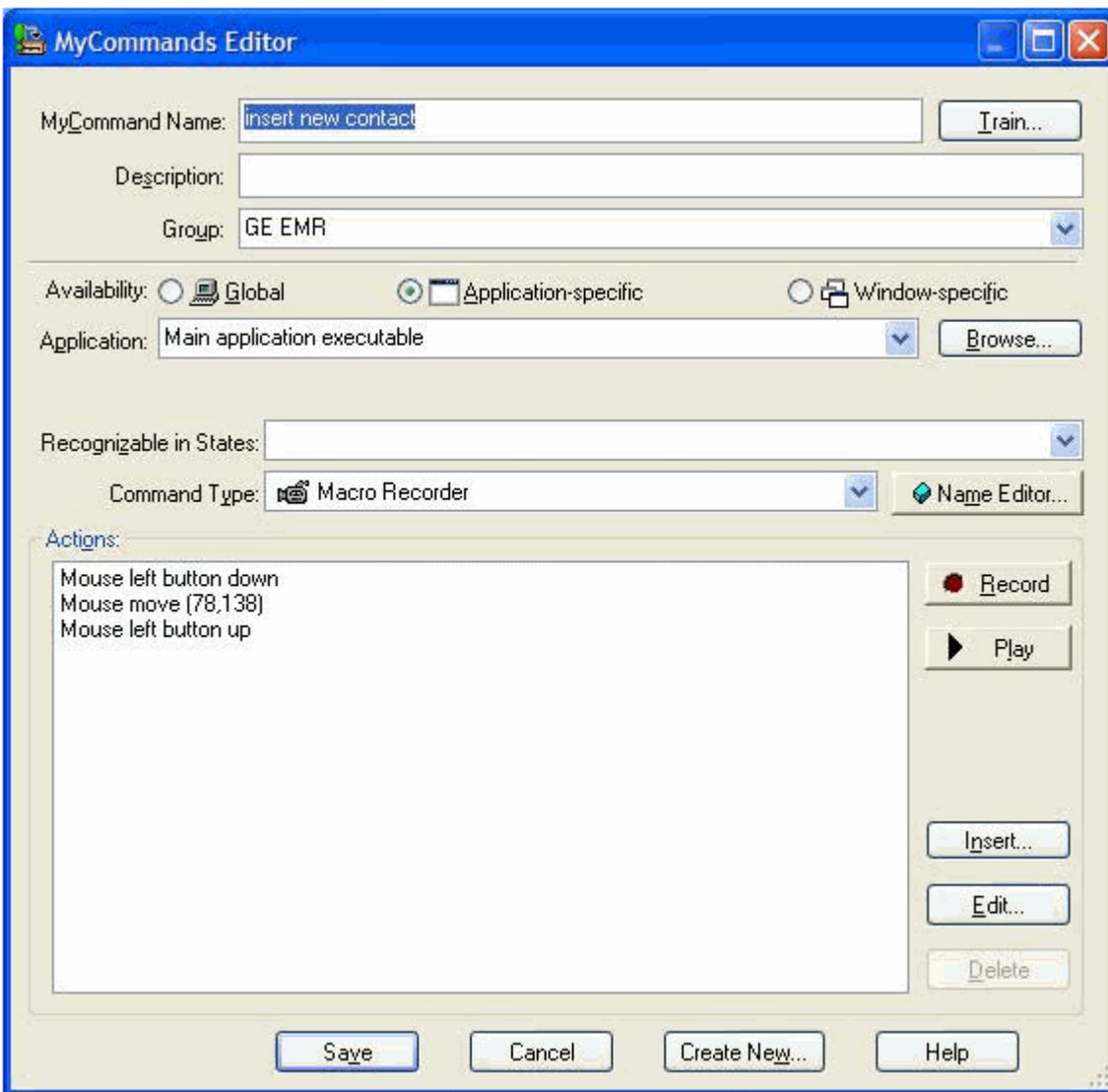

16. **(Optional)** Click the **Play** button to play back the command for testing.
17. When you have finished, click **Save**.

To Use a Macro Recorder Voice Command:

18. Open the application you want to use.
19. With your microphone turned on, say the command name. The actions of the command take place on your screen.